SEA-BAND®

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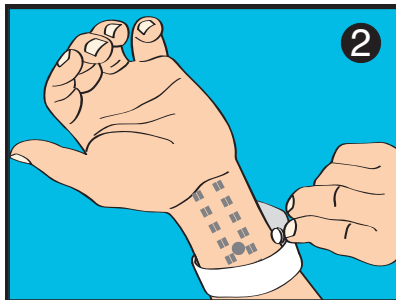
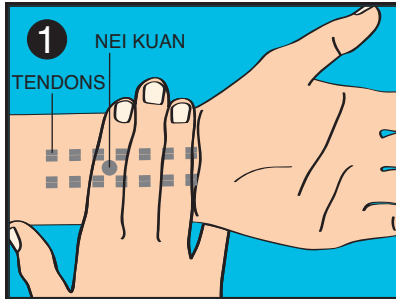
Manufactured and Distributed by:
Sea-Band Ltd, England
Tel: +44 1455 639750
Fax: +44 1455 639751
email: hq@sea-band.com
www.sea-band.com
Sea-Band is a registered trademark

SEA-BAND®
THE NATURAL CHOICE FOR NAUSEA RELIEF

The Sea-Band operates by exerting pressure on the acupressure point (called the Nei-Kuan Point) on each wrist. It takes only a few minutes to become effective.

The Sea-Band provides effective relief from all forms of Travel Sickness without causing drowsiness or other side effects.

It has also been found to be effective in other conditions that produce nausea, such as pregnancy, anaesthesia and chemotherapy.

The Sea-Band may be used repeatedly or worn continuously and can be hand washed using soap or mild detergent. Remove the bands if any swelling occurs.

Sea-Band is suitable for adults and children.

HOW TO POSITION THE SEA-BANDS CORRECTLY

- Place your middle three fingers on the inside of your wrist with the edge of third finger just below the first wrist crease. The Nei-Kuan point is located underneath your first(index) finger in between the two wrist tendons.
- Position the button facing downward over the Nei-Kuan point. Repeat the process for the other wrist.

One band must be worn on each wrist to be effective. Sea-Bands can be put onto your wrist before you start travelling, or after you have started your journey.

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